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
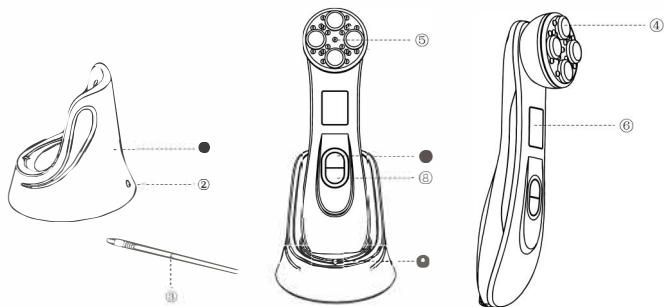
Figure 1
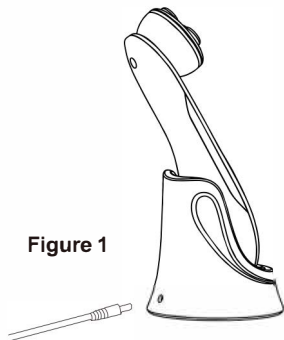
### III OPERATING INSTRUCTIONS

1. Charging – Charge device for 3 hours before first use. To charge, insert the USB cable connector into the device and the other end of the cable into the charging jack and connect to power source.
2. After cleaning skin, apply your preferred serum to skin.
3. Press the Power switch on the device. The red LED indicator will light and the device will be at intensity level 1.

4. Adjust the intensity level according to skin type and personal preference.
5. Adjust intensity by pressing the power button,

- Press once for intensity level 1
- Press twice for intensity level 2
- Press three times for intensity level 3
- Press four times for intensity level 4
- Press five times for intensity level 5
- Press six times to power off the device

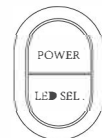
Figure 2

6. Move the device slowly from bottom to top of face and from inside to out. The skin may feel warm during use and this is normal. The device will automatically shut off after 10 minutes of use.

It is recommended that the device be used 3 times per week for 10 minutes each time. It is also recommended that a skin serum be used to keep the skin moist during each treatment with this device.

### IV. LED LIGHT TREATMENT OPERATION

#### Mode 1-Red LED Light Treatment

Red LED light is used to treat wrinkled skin and enhance the regeneration of collagen and skin elasticity. Press the LED selection button once to enter Mode 1.

#### Mode 2-Green LED Light Treatment

Green LED light is used to treat dark spots or discolored skin by working to prohibit the overproduction of melanin. Press the LED selection button twice to enter Mode 2.

#### Mode 3-Blue LED Light Treatment

Blue LED light is used to treat acne and oily skin. Press the LED selection button three times to enter Mode 3.

#### Mode 4-Yellow LED Light Treatment

Yellow LED light is used to brighten the skin. It is suitable for skin with dark spots and dull areas. Press the LED selection button 4 times to enter Mode 4.

#### Mode 5-Pink LED Light Treatment

Pink LED light is used to whiten the skin. Press the LED selection button 5 times to enter Mode 5.

#### Mode 6-Flashing Pink LED Light Treatment

Pink Flashing LED light is used to aid in the permeation of skin care products. Press the LED selection button 6 times to enter Mode 6.

#### Turn off LED Light

Press the LED selection button 7 times to turn off the LED light feature.

Correct operating method  
Metal head fully contact skin

Incorrect operating method  
Metal partial contact skin

## V PARAMETERS AND PACKING INFORMATION

Name: Rejuvenate LED Anti-Aging Device

Standby time:40 days

Input :DC5V 500mA

Use Time:10mins/time

Battery capacity : 750mAh

Unit Weight:220g

Power: Max 10W

Unit Size:165\*45\*46mm

Charging time : 3.5 hours

Box:215\*152\*85mm

## VI SAFETY INSTRUCTIONS

1. Read all safety instructions before using and use correctly.
2. Follow the operating instructions to use this device.
3. Do not pull down or remake the device, please consult the distributor for repairing.
4. Only use the device for its intended purpose as described in this manual.
5. Do not store or use the device in the high temperature and humidity place such as under sunshine or near fire or in the bathroom, etc. To avoid fire or electric shock, Do not use or place it in or near the high temperature appliance.
6. Do not use and place the device in the humid condition, in case of any short circuit, electrical shock or fire.
7. Keep away from children and physical disabled people allow to use, close supervision is necessary when this device is operated by user under 18years of age.
8. Keep out of reach of the children.
9. Pregnant, and people with hypertension, heart disease are not allow to use.
10. The person who is sensitive to the skin diseases excessively is not allowed to use.

## WARNING

1. To avoid malfunction or fire, please stop use in case of unusual hot.
2. Please do not drop or clash the device.
3. Please apply serum constantly during use.
4. Please do not use it if the power cord is damaged.
5. Below conditions are not allow to use:  
Heart disease, fever, communicable disease, cancer, open cuts, acute disease, Hypotension, hypertension, install pacemaker, artificial heart or lung, during medical treatment, skin disease, sunburned, over sensitive people, anapetia.